

An SP Industries Company

Micro-Vial Homogenizer System

“Pellet Pestles®”

**Introducing the State-of-the-Art System
for reconstituting pelleted materials or cells in microfuge tubes.**

The Wilmad-LabGlass Micro-Vial Homogenizer System features an advantaged microfuge tube that is specially matched to a micro pestle which provides the tightest tolerances between the vial wall and the pestle. The result: the best homogenization of precipitates or cells in the industry.

pestle



0.5 mL Tube



1.5 mL Tube



Key Features

- Vials and Pestles made of animal and phalate/phalate-ester free polypropylene with lot control and material traceability.
- Certified DNase, RNase and Pyrogen-free.
- Colorant-free to eliminate potential leaching of heavy metals into samples.
- Vials have 25% thicker bottoms to ensure integrity during use.
- Unique ribbed shaft provides precise finger control of homogenization for shear sensitive samples.
- Matched vials and pestles available in 0.5 ml & 1.5 ml sizes.
- Available individually-packaged or in a unique bulk-pack style.
- Optional Battery-powered Cordless Motor available to process larger numbers of samples or routine applications.



Phone

Toll Free: **1-800-220-5171**
Intl.: **856-691-3200**




Fax

Toll Free: **1-800-220-1081**
Intl.: **856-691-6206**



Web

www.wilmad-labglass.com
cs@wilmad-labglass.com

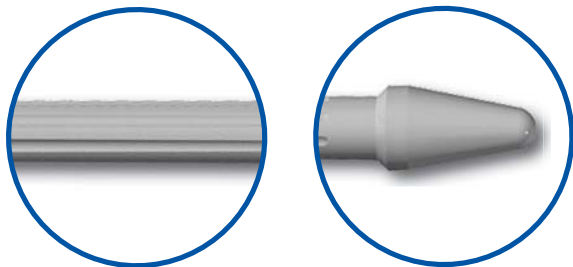

Micro-Pestle Features

- Ribbed shaft provides improved finger control without motor.
- Allows air to escape for easier insertion into motor adapter.
- Matched contact angle of tip to vial wall to ensure the best reconstitution of materials.
- Offset parting line in mold design eliminates particulates often generated during use of other homogenizers in the marketplace.



Micro-Centrifuge Vial Features

- 0.5 ml & 1.5 ml Sizes
- Positive seal "Snap-Cap" Design
- Non-graduated, with frosted marking spot for easy sample identification
- Flat frosted cap surface for convenient labeling
- Pierceable Top



Packaging Features

- Style 1) Individually-wrapped vials and pestles over packed in resealable bags of 100 each.
- Style 2) Unique bulk-packaged products supplied in a resealable bag with an additional heat-seal above the standard closure that indicates if the package has been opened.

Individual Pack, Autoclavable

Product Number	Description	Case Qty.	Price
BP-7001-005	0.5 mL Pestle	100	\$95.00
BP-7002-005	0.5 mL Tube	100	\$60.00
BP-7000-005	0.5 mL Pestle and Tube	100	\$105.00
BP-7001-015	1.5 mL Pestle	100	\$95.00
BP-7002-015	1.5 mL Tube	100	\$60.00
BP-7000-015	1.5 mL Pestle and Tube	100	\$105.00

Bulk Pack, Autoclavable

Product Number	Description	Case Qty.	Price
BP-7011-005	0.5 mL Pestle	100	\$74.00
BP-7012-005	0.5 mL Tube	100	\$42.00
BP-7010-005	0.5 mL Pestle and Tube	100	\$89.00
BP-7011-015	1.5 mL Pestle	100	\$74.00
BP-7012-015	1.5 mL Tube	100	\$42.00
BP-7010-015	1.5 mL Pestle and Tube	100	\$89.00



Product Number	Included Items	Price
BP-7005-000	<ul style="list-style-type: none"> ✓ Cordless Motor ✓ 2 AA Batteries ✓ Pestle Adapter 	\$78.00



Toll Free: **1-800-220-5171**
Intl.: **856-691-3200**



Toll Free: **1-800-220-1081**
Intl.: **856-691-6206**



www.wilmad-labglass.com
cs@wilmad-labglass.com

Pellet Pestle® is a registered trademark of Kontes (Kimble-Chase)